


Metabolic Cage For Large Animal (Suitable for Rabbit & Guinea Pig)

FOR SEPARATION OF FAECES AND URINE

Orchid's metabolic cages for large rodents are equipped with system for separation of faeces and urine. Each system is divided in 3 parts- cage body, separation hopper and feed- water supplying unit. Separation hopper is fitted with S.S. mesh for separation of urine and faeces. Faeces are collected above the mesh while urine goes down through the tube. Design and finish are such that the animal has a maximum comfort. Maintenance is efficient and fast. Every unit comes with cage floor, food hopper, water bottle, urine collection system.



**Designed To Offer Better
Separation Efficiency Of
Urine And Feces Assuring
Maximum Purity Of Samples**

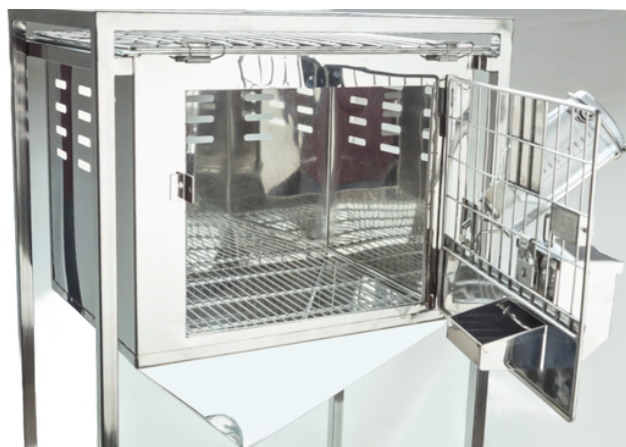

COMPONENTS INCLUDED:

- Made from complete AISI S.S. 304 quality material
- Overall Dimensions along with stand: 535X390X825mm (LXWXH)
- Dimensions of cage body with funnel: 480X390X610mm
- Cage Stand is made from 20X20X15mm S.S. 304 square pipe
- Cage body and funnel is made from S.S 304 sheet of 0.8mm thickness
- Top grill is made from 2.5mm and 5mm S.S. 304 bar has gap of 25X60mm
- Door has dimensions of 330X280mm and is made from 2.5mm and 5mm S.S. 304 bar has gap of 25X45mm
- Floor grill is made from 2.5mm and 5mm S.S. 304 bar has gap of 11X80mm

ORDERING INFORMATION:

Product Name	Model	Accessories Supplied
Metabolic Cage For Large Animal	DUC-02	<ul style="list-style-type: none">• Cage body with funnel• Floor grill, stand• water bottle• Bottle holder• Waster water collection tube• Food hopper• S.S separation plate• Urine collection bottle of 500ml



Orchid Scientific & Innovative India Pvt. Ltd.

📍 B-59, M.I.D.C., Ambad, Nashik - 422010, India.

☎ +91253-2387600, 2972525

✉ office@orchidscientific.com, exports@orchidscientific.com

🌐 www.orchidscientific.com



www.orchidscientific.com