

Metabolic Cages For Mouse


METABOLIC CAGES FOR MOUSE

(ECONOMY MODEL- DIURESIS CAGE WITH STAND):

MODEL: DUC-1-M

Especially useful for urine collection experiments. The cage is provided with cage body, perforated floor, Funnel, drinking water bottle of capacity 15 ml, waste water collection tube and urine collection tube. Cage body is fitted with tag holder for identification purpose. Urine and faeces falls on the S.S. Funnel and one separation plate of S.S. is fixed inside funnel which allows passage of only urine and not faeces, leading to collection of faeces above the separation plate and urine at the bottom in the P.P. tube.

Floor space: 94 sq. cm.

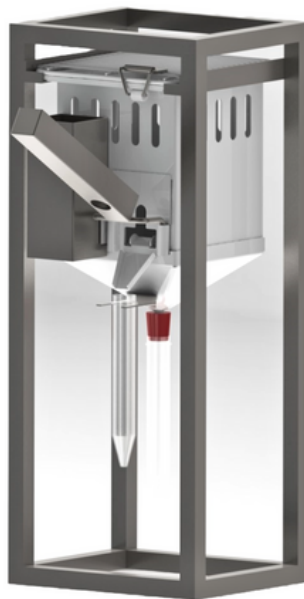

COMPONENTS INCLUDED:

- S.S. Cage body with funnel
- S.S. perforated floor
- S.S. separation plate for separation of urine and faeces
- Water bottle with rubber cork of capacity 15ml: 01
- Food hopper: 01
- Water bottle holder: 01
- Waste water collection P.P. tube: 01
- Urine collection tube P.P. with rubber cork: 01
- S.S. Stand for holding single cage: 01

Overall Approx. Dimensions: 100 x 100 x 120mm (L X W X H)



**Designed To Collect Urine
And Waste Water. Especially
Useful For Urine Collection
Experiments.**

**METABOLIC CAGES FOR MOUSE
(ECONOMY MODEL- DIURESIS CAGE WITHOUT STAND, SUITABLE FOR RACK):**

MODEL: DUC-1-M-WS

Especially useful for urine collection experiments The cage is provided with a cage body, perforated floor, Funnel, drinking water bottle of capacity 15ml, waste water collection tube and urine collection tube Cage body is fitted with tag holder for identification purpose. Urine and faeces falls on the S.S. Funnel and one separation plate of S.S. is fixed inside funnel which allows passage of only urine and not faeces, leading to collection of faeces above the separation plate and urine at the bottom in the P.P. tube.

Floor space: 94 sq. cm.

COMPONENTS INCLUDED:

- S.S. Cage body with funnel
- S.S. perforated floor
- S.S. separation plate for separation of urine and faeces
- Water bottle with rubber cork of capacity 15 ml
- Food hopper: 01
- Water bottle holder: 01
- Waste water collection P.P. tube: 01
- Urine collection tube P.P. with rubber cork: 01

RACK FOR METABOLIC CAGES (SUITABLE FOR DUC-1-M-WS):

Made of S.S. 304 material and suitable to hold 6-8 metabolic cages for small rodents (economy model) DUC-1-M-WS Fitted with castor wheels out of which two are with break and two are without break



ORDERING INFORMATION:

Model	Outer dimensions approx. with wheel (mm)	Number of shelves	Number of cages/shelf	Total cages/trolley
DUC-M-SS-15	840 x 450 x 1527	3	5	15
DUC-M-DS-30	840 x 450 x 1527	3	5 (Double side)	30



Orchid Scientific & Innovative India Pvt. Ltd.

📍 B-59, M.I.D.C., Ambad, Nashik - 422010, India.

☎ +91253-2387600, 2972525

✉ office@orchidscientific.com, exports@orchidscientific.com

🌐 www.orchidscientific.com



www.orchidscientific.com